

28U25 WalkOn Reaction junior

It's so dynamic!



WalkOn Reaction junior

Dynamic ankle foot orthosis
made of carbon for children

The WalkOn Reaction junior is the new member of the WalkOn line. Based on the WalkOn Reaction plus, it features highly dynamic properties. The ankle foot orthosis uses the ground reaction force to significantly improve the user's gait pattern. It is geared towards children, who require greater support than that provided by traditional dorsiflexion assist orthoses.



Benefits at a glance

- Highly dynamic properties (use of ground reaction force)
- Carbon fibre prepreg material for high energy return
- High wearer comfort due to low weight and open heel and the anatomically designed frontal support element
- Shin pad made of temperature - regulating material for improved skin hygiene and wearer comfort
- Combination with additional dynamic ankle foot orthoses (DAFOs) possible
- All fabric parts can be washed and replaced
- Sole can be trimmed to size
- Can be worn with sturdy standard shoes

Indications

- Cerebral palsy
- For dorsiflexor weakness, without or with mild to moderate spasticity
- Hemiparesis, diparesis

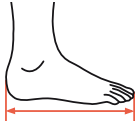
Effects

- Supports lifting of the foot
- Promotes dynamic rollover
- Can help increase the stride length by improving muscle dynamics
- Helps promote exterior rotation of the foot at heel strike
- Helps promote a physiological gait pattern
- May reduce the risk of falling





The orthosis design is adapted to the anatomical conditions of children and for this reason differs in its various sizes. Because it is comfortable to wear, the WalkOn Reaction junior enjoys high acceptance amongst children.

Art. no. 28U25			
Order no.: Article number=side & shoe size Example for ordering: 28U25=L24-27			
	Side	Shoe size	Overall height (cm)
	L/R	24-27	20.5
	L/R	27-30	24.5
	L/R	30-33	27
	L/R	33-36	29.4

A FIT KIT (28T5) is available for the WalkOn Reaction junior, as it is for the other products in the WalkOn line.

© (2016) Tabaluga Enterprises GmbH.

TABALUGA is based on an idea from
Peter Maffay, Gregor Rottschalk and
Rolf Zuckowski. Images by Helme Heine.

www.tabaluga-enterprises.de

Ottobock Australia Pty Ltd
1.01, 62 Norwest Blvd · Baulkham Hills NSW 2153
T +61 2 8818 2800 · F +61 2 8814 4500
healthcare@ottobock.com.au · www.ottobock.com.au